Administration Training for School Staff

Plainville School District

Massachusetts Department of Public Health





Purpose

To provide unlicensed school personnel with the basic knowledge and skills to administer epinephrine by auto-injector in a lifethreatening situation.



DPH Regulations

- School Registration by nurse with the Department of Public Health(Renew every 2 years)
- Initial training for competency
- Training review 2 times yearly



Outline

Watch Food Allergy training module on the web:

Food Allergies and Anaphylaxis in School: What School Staff Need to Know

- Continue with the power point:
 - Watch 3 video clips on administration of each type of epipen
 - Review School Emergency Plans
- Competency
 - Competency test
 - Make appointment with nurse to demonstrate administration



Stop here and watch the training video....click on link below:

Epinephrine Training Video

Once have watched all modules then continue on...



Administration of Epipen

 3 types of epinephrine auto-injectors available



1. Epipen

2. Adrenaclick

3. Auvi-Q

EpiPen Auto-injector Yellow case-Adult dose **Green** case-Junior dose



- Pop open the yellow or green cap off the carrying case and remove EpiPen from its storage tube
- 2. Grasp unit with orange tip pointing downward
- 3. Form fist around the unit (orange tip down).
- 4. With your other hand, pull off the blue safety release cap.
- 5. Swing and jab firmly into outer thigh until it clicks so that unit is at 90 degree angle with thigh. (can go through clothing)
- 6 Hold firmly against thigh for 10 seconds
- 7. Remove unit from thigh and massage injection site for 10 seconds
- Call 911 and seek immediate medical attention. Symptoms may reoccur.



Click on link below to watch a demonstration video on how to administer Epipen

(halfway down web page)

Epipen Demonstration



Adrenaclick Auto-injector






- 1. Pull off grey cap labeled "1": you will now see a red tip. put hand or fingers over the red tip.
- 2. Pull of grey end cap labeled "2".
- 3. Place the red tip against the middle of the outer side of the thigh (upper leg). It can go through clothes.
- 3. Press down hard until auto-injector fires and the needle enters the thigh. Hold in place while slowly counting to 10.
- 4. Remove Adrenaclick from thigh and check the red tip. If the needle is exposed, the dose was given.
- 5. Immediately call 911. Symptoms may reoccur.



 Click the link below to see a video on how to administer Adrenaclick:

Adrenaclick Demonstation

(Choose caregiver administration video)

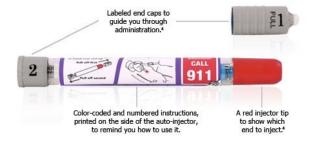

Auvi-Q auto-injector

- 1st auto-injector that has a voice that talks you through each step of the injection process
- Counts down when you are done, beeps and has a red flashing light to tell you injection was successful
- Compact Design fits in pocket
- 2 doses
 - * Child dose 0.15mg (33-66 pounds)
 Blue packaging
 - * Adult Dose 0.3mg (over 66 pounds)

 Orange packaging
- Needle injects and retracts automatically





 Click the link below to watch a video on how to administer Auvi-Q

Auvi-Q Demonstration

(video may have a short delay to load)





School Emergency Plan

- Notify Nurse and main office of episode immediately!!
- If Nurse unavailable-call main office and they will institute the emergency plan.
- Red Emergency Care Plan Book is in the Main office by Enrollment info-has copies of all student emergency plans



Emergency Plan Steps for Trained Staff

If the nurse is unavailable you will be notified of need for Epinephrine administration:

- 1. Obtain the Epinephrine from the nurse's office, cafeteria or on student.
- 2. Administer the epinephrine
- 3. Notify **911** immediately/designate someone to do this-student must be transported to the ER even if symptoms improve after treatment.
- 4. Notify Principal and Parent-emergency contact info on back of emergency plan.
- 5. Department of Public Health to be notified by Nurse



Competency Award

Post-test on Web. Click on the link below:

Competency Test

 Call nurse to set up time to demonstrate the administration of the Epinephrine auto-injector

CONGRATULATIONS!